

Pigtail Adapter


Documentation v0

John Hermiz

11/7/16

Overview

- This document describes how to use, take apart and assemble the pigtail adapter

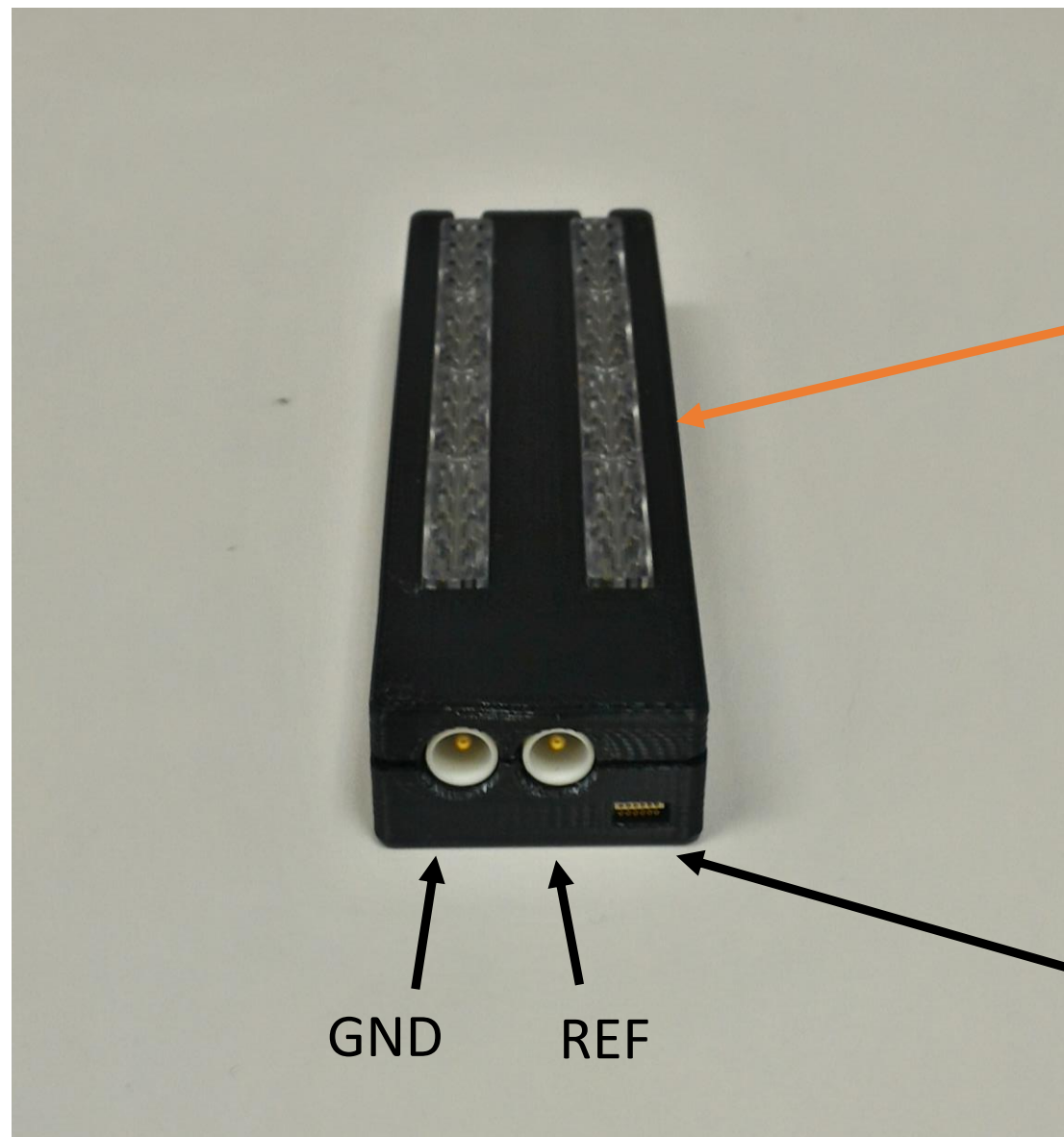


Pigtail Connector

GND

REF

Intan Connector

How to Setup

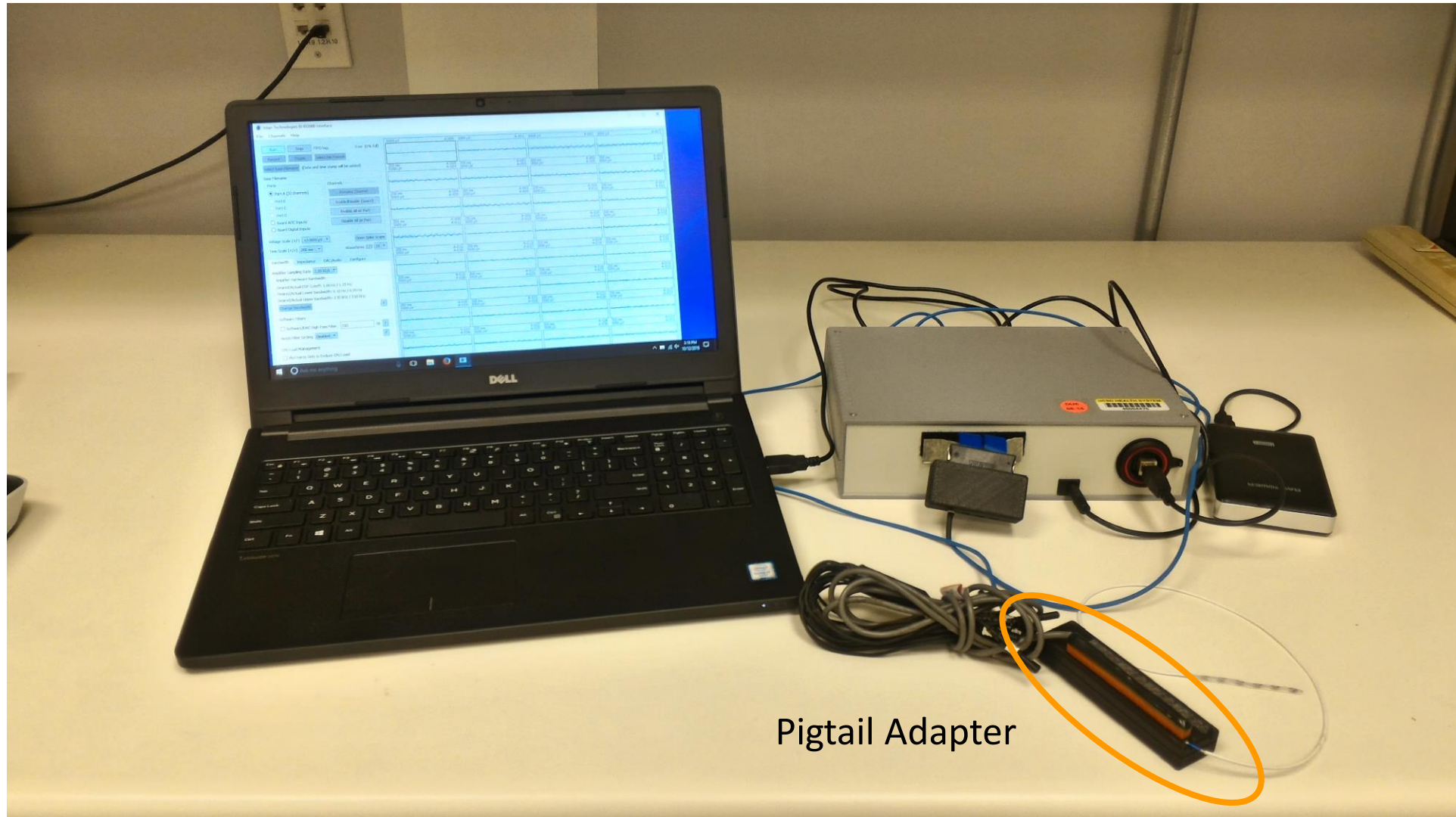


Pigtail Adapter

This is how the final setup looks

Connecting

1. Connect blue Intan cable (for research hardware)
2. Connect Gnd and Ref
3. Slip pigtail into connector
 1. Ensure it's properly inserted!
4. Connect clinical cable (for clinical hardware)

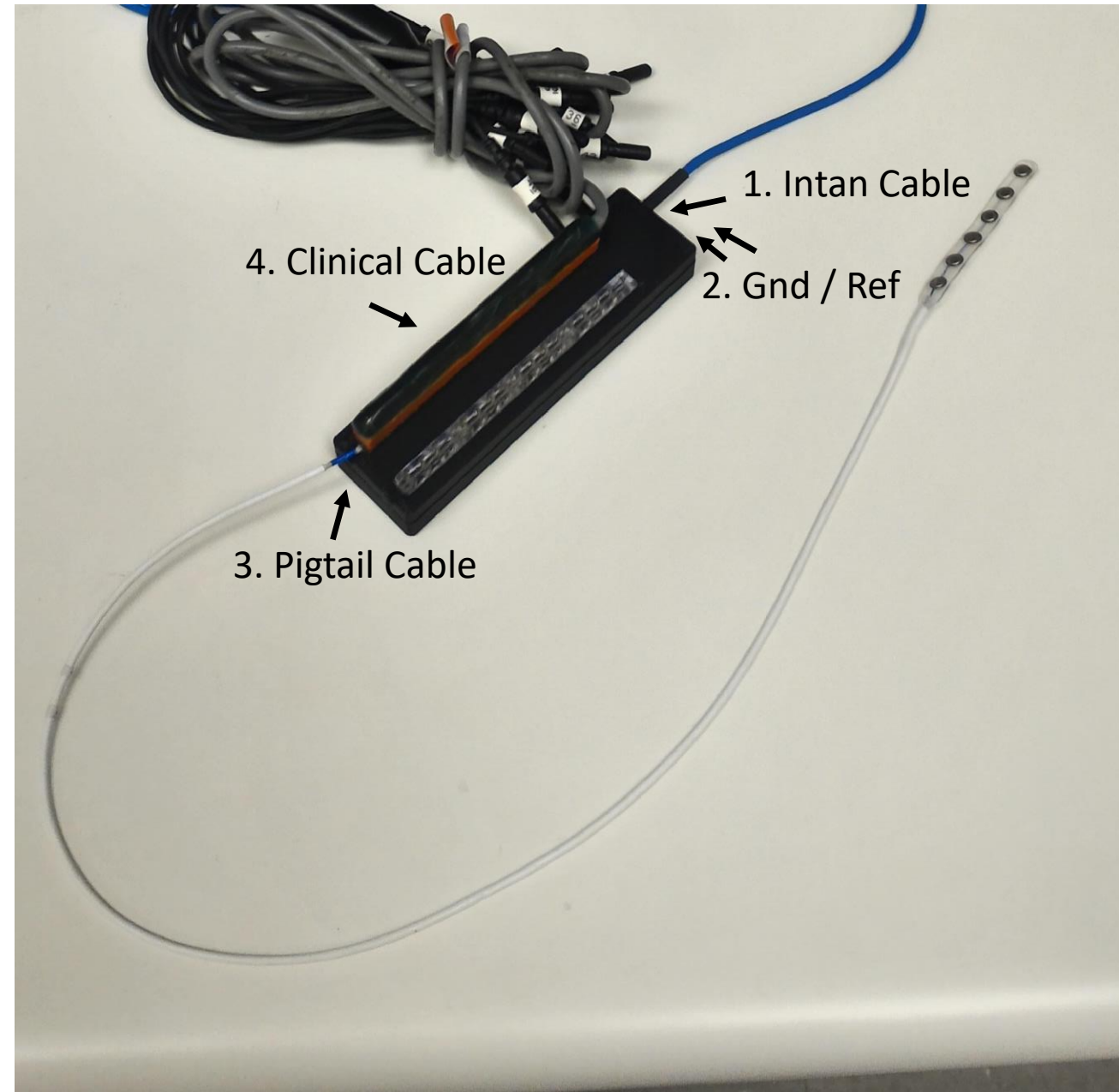


Dis-assemble

WATCH DISASSEMBLE VIDEO

Written Instructions

1. Remove all connector
2. Use Philips to remove 2 screws from back
3. Pull apart back and front of case
4. Carefully pull away electronics
5. Remove pigtail connectors by pushing them inwards

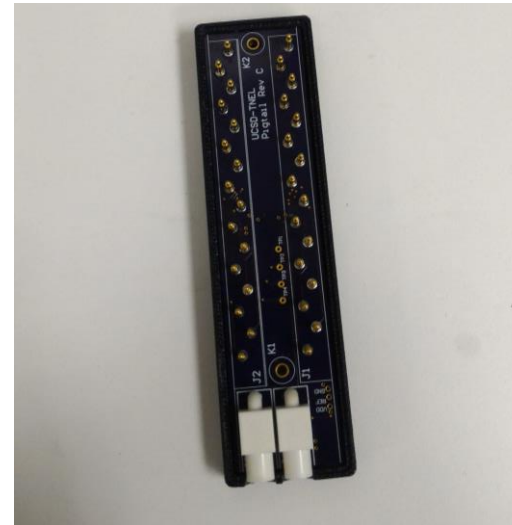

Electronics board



Pigtail Connectors can be removed for cleaning or replacement

Assemble

WATCH ASSEMBLE VIDEO

Do disassemble instructions in reverse



Push pigtail connector back into place